



29th Annual Meeting of the
American Society for Gravitational and Space Research (ASGSR)

November 3rd – 8th, 2013
Orlando, Florida, USA


MICROGRAVITY SCIENCE GLOVEBOX (MSG) SPACE SCIENCE'S PAST, PRESENT, AND FUTURE ON THE INTERNATIONAL SPACE STATION (ISS)

Reggie A. Spivey

Teledyne Brown Engineering
Huntsville, Alabama

Lee P. Jordan

NASA Marshall Space Flight Center,
Huntsville, Alabama

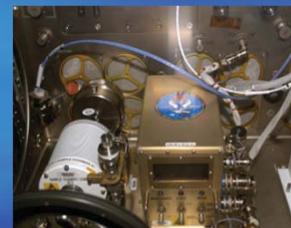



Microgravity Science Glovebox



Agenda

- Introduction
- Payload Interfaces and Resources Provided by MSG
- Overview of the Research Accomplished in the MSG Facility to Date
- MSG Operations Planned for 2014
- Life Science Ancillary Hardware (LSAH) Upgrades
- Video Upgrade Equipment (VUE)
- Conclusion





Microgravity Science Glovebox



Introduction

- The Microgravity Science Glovebox (MSG) is a double rack facility designed for microgravity investigation handling aboard the International Space Station (ISS).
- The unique design of the facility allows it to accommodate science and technology investigations in a “workbench” type environment
- MSG facility provides an enclosed working area for investigation manipulation and observation in the ISS. Provides two levels of containment via physical barrier, negative pressure, and air filtration .
- The MSG team and facilities provide quick access to space for exploratory and National Lab type investigations to gain an understanding of the role of gravity in the physics associated research areas.





Microgravity Science Glovebox

Microgravity Science Glovebox (MSG)



Marshall Space Flight Center

MSG Facility Hardware Overview

Removable Side Ports

16" diameter on both Left and Right sides for setting up hardware in Work Volume


Glove Ports

Four identical glove ports are located on the left and right side loading ports and the front window

DC Power Switching And Circuit Breakers

Stowage Drawers

Video System Drawer

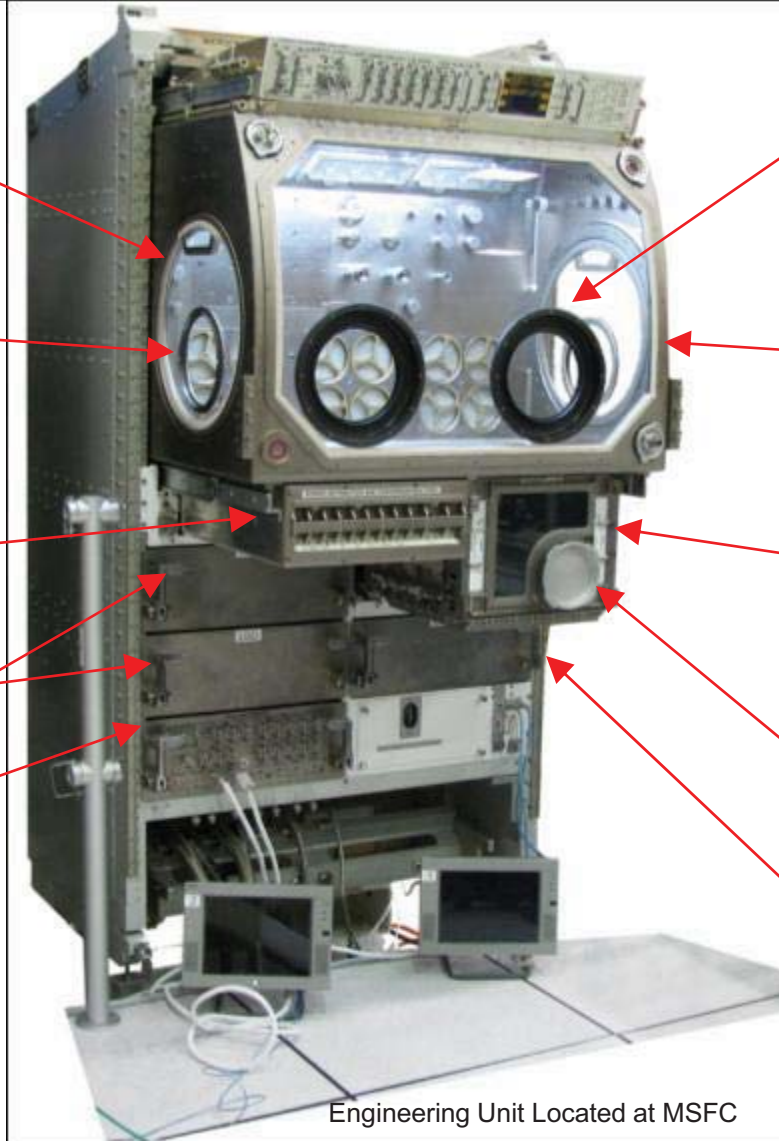

Front Window Glove Ports

Four 6" diameter glove ports can be fitted with any of three different sized gloves or blanks

Core Facility

Retractable Core Facility includes the Work Volume, Airlock, Power Distribution & Switching Box, and the Command and Monitoring Panel

Airlock

Provides a "Pass Through" for hardware to enter the Work Volume without breaking Containment. The lid of the Air Lock opens up into the floor of the Work Volume

Airlock Glove Port with Blank

A Single 4" diameter glove port can also be fitted with any of three different sized gloves or a blank

Stowage Drawers

Engineering Unit Located at MSFC



Microgravity Science Glovebox

Microgravity Science Glovebox (MSG)

Current MSG-Provided Payload Interfaces/Resources



Marshall Space Flight Center

-
- **Work Volume(WV) - Volume**
 - 0.255 m³ = 255 liters
 - **Work Volume - Dimensions**
 - 906mm wide x 637mm high
 - 500mm deep (at the floor)
 - 385mm deep (at the top)
 - **Maximum size of single piece of equipment in WV (via side access ports)**
 - 406mm diameter
 - **Payload Attachment**
 - M6 threaded fasteners in floor, ceiling, & sides
 - **Power available to investigation**
 - +28V DC at useable 7 amps
 - +12V DC at useable 2 amps
 - -12V DC at useable 2 amps
 - +5V DC at useable 4 amps
 - +120V DC at useable 8.3 amps
 - **Maximum heat dissipation**
 - 1000W Total
 - 800W from coldplate
 - 200W from air flow
 - **General illumination**
 - 1000 lux @ 200mm above WV floor
 - **Video**
 - 4 color Hitachi HV-C20 cameras
 - 2 Sony DSRV10 Digital Recorders
 - 2 Sony GV-A500 Analog 8mm Recorders
 - **Data handling connections**
 - Two RS422-to-MSG for investigations
 - One MIL-BUS-1553B-to-MSG for communication via MLC
 - Ethernet LAN 1 and LAN 2 (in US LAB)
 - MSG Laptop Computer (MLC) – IBM T61P
 - **Filtration**
 - 12 HEPA/charcoal/catalyst WV filters
 - **1 HEPA/charcoal/catalyst Airlock filter**
 - **Up to Two Levels of Containment**
 - Physical barrier of MSG structures, gloves, etc.
 - Negative pressure generated by MSG fans.
 - **Other resources available**
 - Gaseous Nitrogen
 - Vacuum (VRS & VES)




Microgravity Science Glovebox

Microgravity Science Glovebox (MSG)



MSG Flight Unit Cumulative Hours of Operation

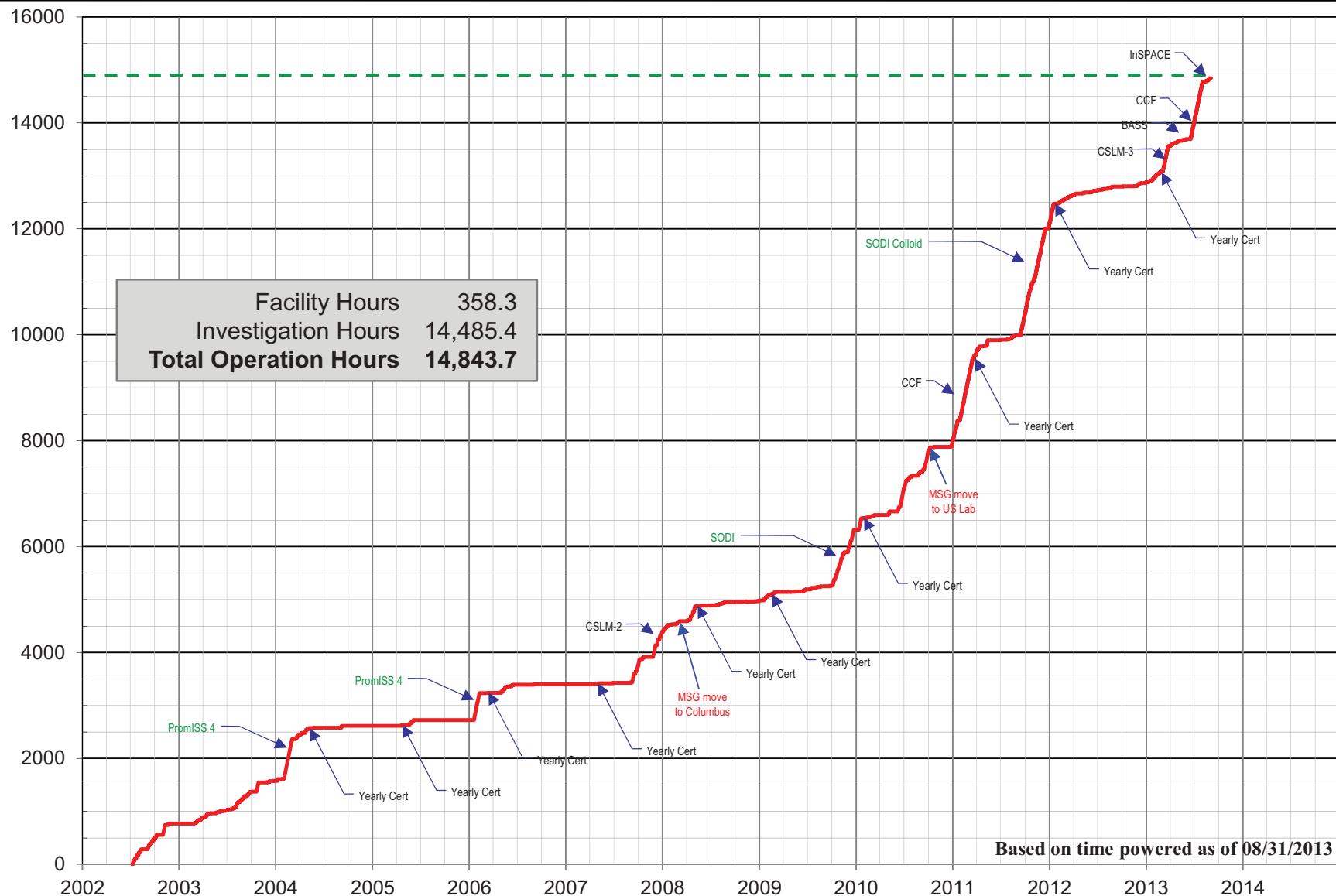





Microgravity Science Glovebox

Microgravity Science Glovebox (MSG)

International Utilization of the MSG Facility

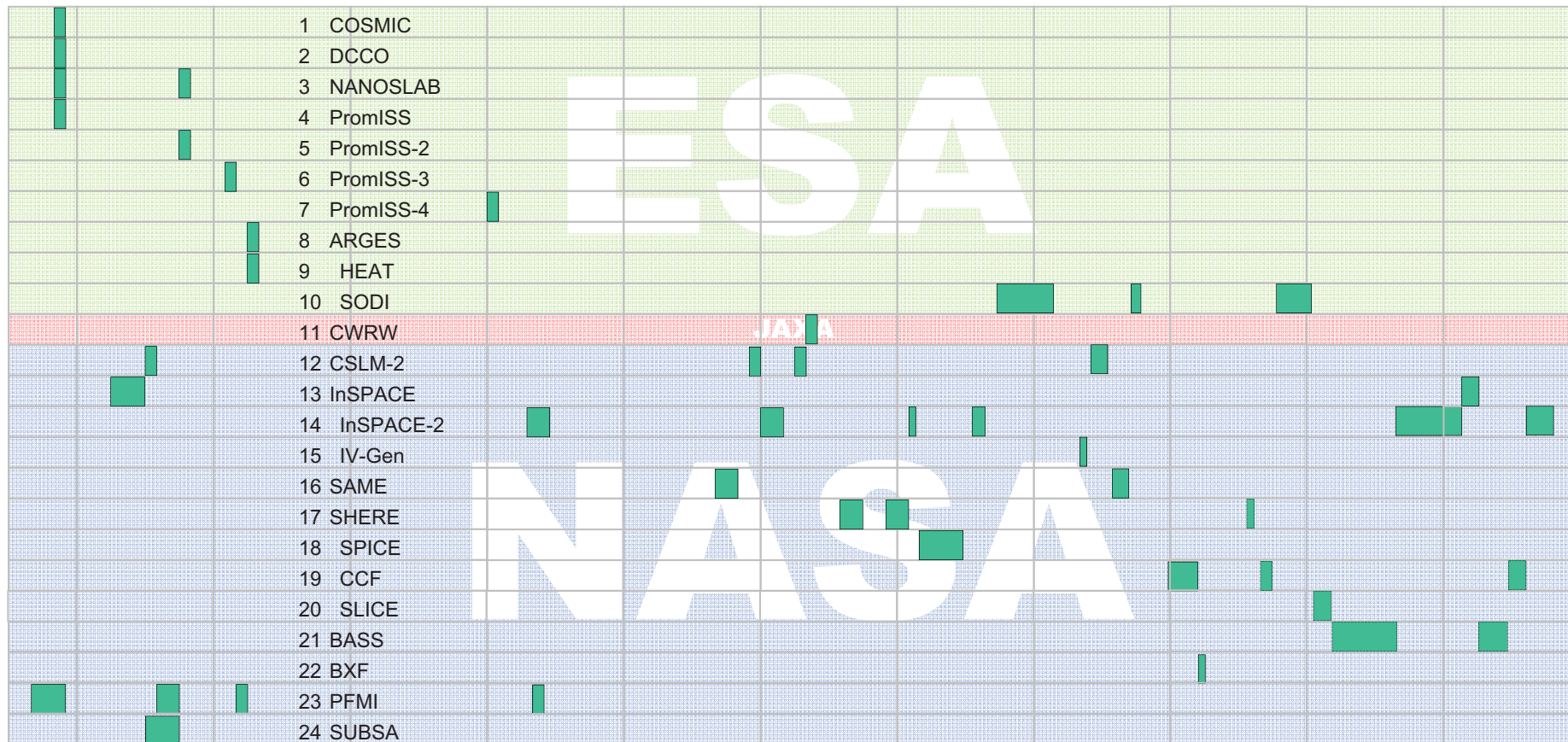

as of 08/31/2013

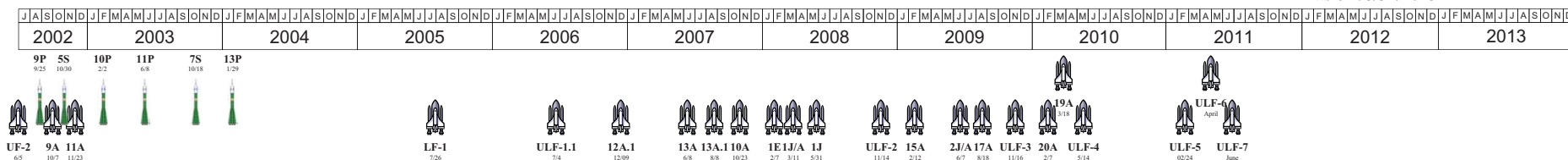



Microgravity Science Glovebox

Microgravity Science Glovebox (MSG)

MSG Investigations



Payload Name & Acronym	Sponsoring Organization	Type of Investigation
Combustion Synthesis under Microgravity Conditions (COSMIC)	ESA	Combustion
Microgravity Experiment for the Measurement of Diffusion Coefficients in Crude Oil (DCCO)	ESA	Diffusion
NANOSLAB	ESA	Zeolite Crystal Growth
Protein Microscope for the International Space Station (PromISS-1,2,3, & 4)	ESA	Protein Crystal Growth
ARGES	ESA	Light Bulb Technology
HEAT	ESA	Heat Pipe Technology
Selectable Optical Diagnostics Instrument (SODI)	ESA	Diffusion and Soret Phenomena
Cell Wall/Resist Wall (CWRW)	JAXA	Plant Growth
Coarsening in Solid Liquid Mixtures-2 (CSLM-2)	NASA	Material Science
Investigating the Structure of Paramagnetic Aggregates from Colloidal Emulsions (InSPACE-1,2, & 3)	NASA	Magnetorheological (MR) Fluids
IntraVenous Fluids GENeration and mixing (IV-Gen)	NASA	Human Health
Smoke Aerosol Measurement Experiment (SAME)	NASA	Spacecraft Smoke Detection
Shear History Extensional Rheology Experiment (SHERE)	NASA	Polymer
Smoke Point Coflow Experiment (SPICE)	NASA	Combustion
Critical Velocities in Open Capillary Channels (CCF)	NASA	Fluids
Structure and Liftoff in Combustion Experiment (SLICE)	NASA	Combustion
Burning and Suppression of Solids (BASS)	NASA	Combustion
Boiling eXperiment Facility (BXF)	NASA	Heat Transfer
Pore Formation and Mobility Investigation (PFMI)	NASA	Material Science
Solidification Using a Baffle in Sealed Ampoules (SUBSA)	NASA	Material Science
Rodent Research	NASA	Life Science
3D Printer	NASA	Technology Demonstration
Bioculture Systems	NASA	Life Science
Observation and Analysis of Smectic Islands in Space (OASIS)	NASA	Material Science
Zero Boil-Off Tank (Z-BOT)	NASA	Heat Transfer
Packed Bed Reactor Experiment (PBRD)	NASA	Physical Science
Transparent Alloys	ESA	Material Science




Microgravity Science Glovebox (MSG)



MSG Operations Planned for 2013-2014

Microgravity Science Glovebox

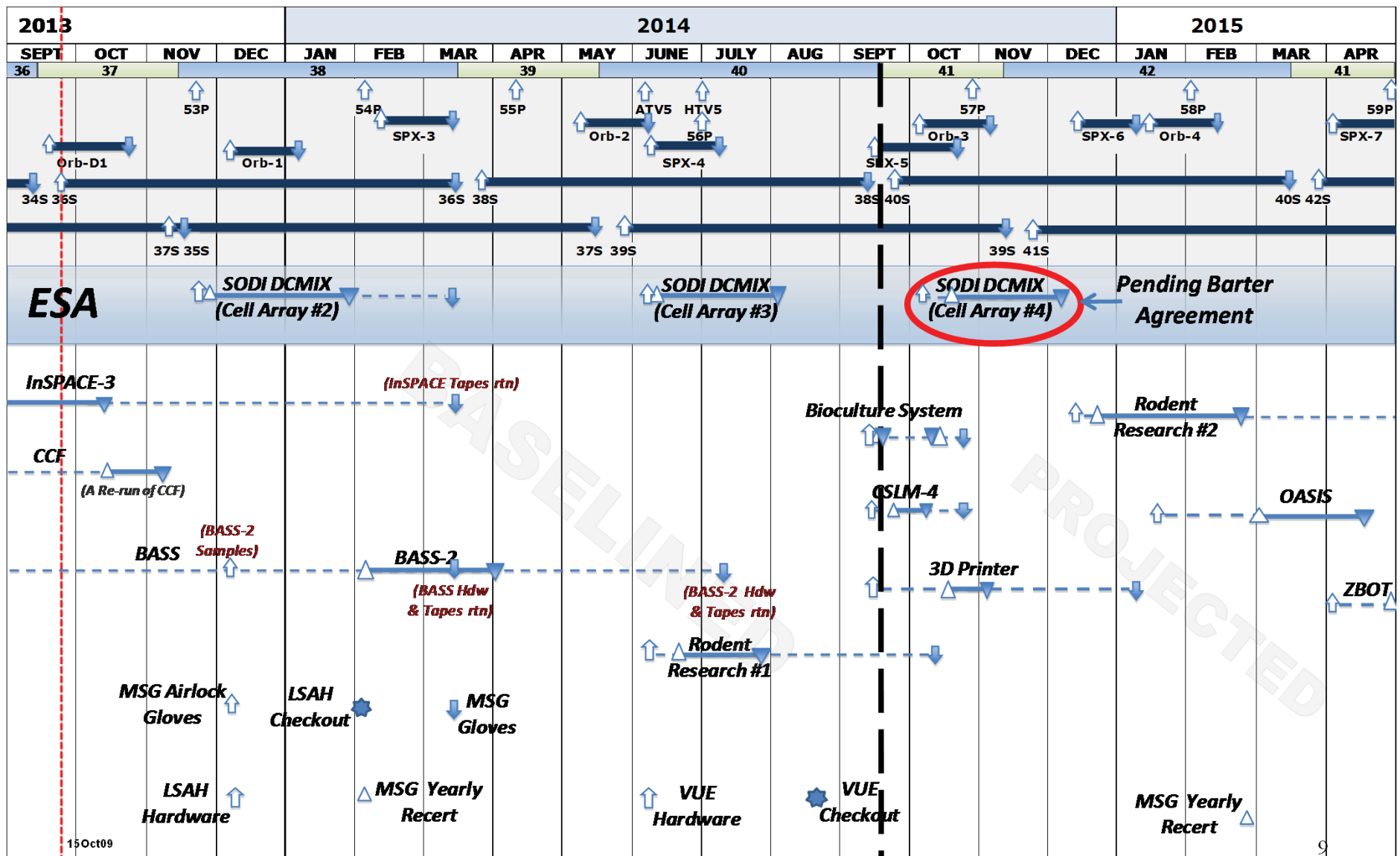



Microgravity Science Glovebox

Microgravity Science Glovebox (MSG)



Marshall Space Flight Center

Life Science Ancillary Hardware (LSAH) Upgrades Available in 2014



Microgravity Science Glovebox

Microgravity Science Glovebox (MSG)



MSG LSAH Upgrades


- Materials utilized by Life Science/Biological Research payloads will require additional capabilities for handling and clean up:
 - Filtration System: a capability added to the existing MSG Work Volume air circulation system that scrubs typical life science biological and chemical contaminants from the MSG Work Volume air.
 - Decontamination System: a capability to reduce released biological contaminants (Bio Safety Levels (BSL) 1 and 2) to levels safe for crew exposure and a capability to remove released contaminants from surfaces within the Work Volume.
 - Exchangeable Glove System this is more suited for various life science activities.



MSG Life Science Filters



Decontamination System

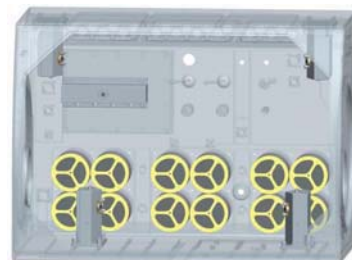

Glove & Gauntlet Configuration



Iris & Gauntlet w/Disposable Glove






Microgravity Science Glovebox

Microgravity Science Glovebox (MSG)



Biological Filters

- MSG's Air Handling Unit creates negative pressure in the Work Volume to provide one means of containment
 - Filter banks trap contaminants when air passes once through the filters
 - Current filter components trap typical material-science and combustions contaminants
- New filters will be added to the existing MSG filters
- New filters will trap typical life/biological science contaminant/materials
 - Such as preservatives, fixatives, and other byproducts

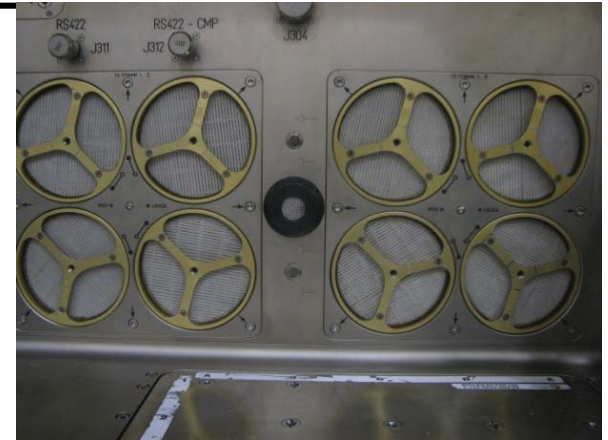

MSG Life Science Filters



**Sundstrom SR 299-2 ABEK1HgP3R
Combination Filter**



In MSG's current design, each of the thirteen front filters is easily exchangeable on orbit by the crew.






Microgravity Science Glovebox

Microgravity Science Glovebox (MSG)




Decontamination System

- New Decontamination Capability within MSG Work Volume
 - Decontaminate before experiment to prevent contamination of biological samples
 - Decontaminate after experiment to disinfect any released biological materials
- Ground-based labs typically use UV Light or Ozone

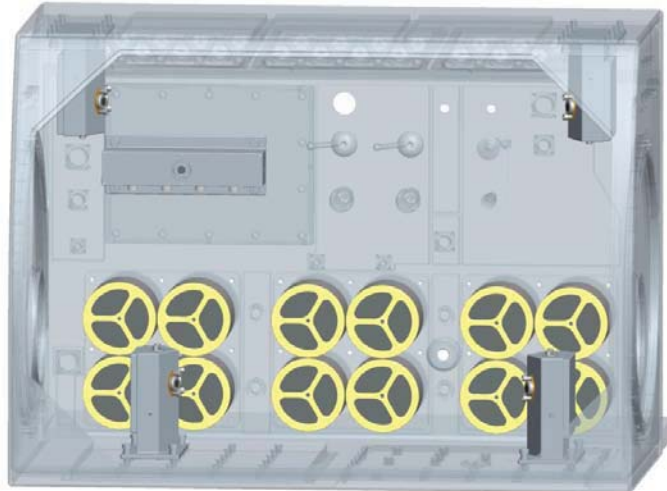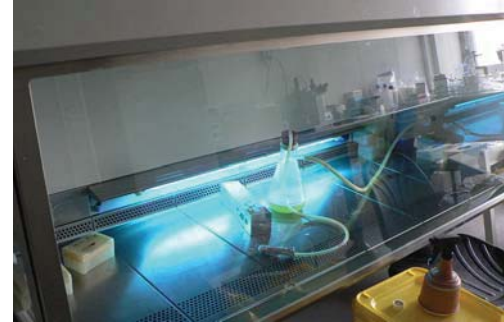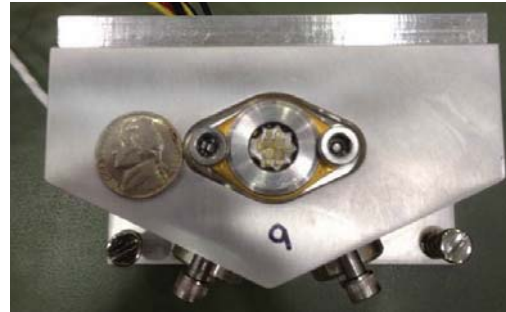

MSG Decontamination System

Ultraviolet germicidal irradiation is a sterilization method that uses ultraviolet light at sufficiently short wavelength to break down microorganisms. It is used in a variety of applications, such as food, air and water purification.

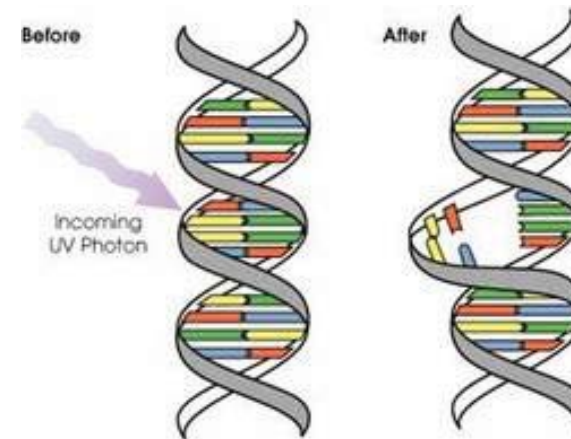



Microgravity Science Glovebox

Microgravity Science Glovebox (MSG)

Decontamination System



Marshall Space Flight Center

List of Microorganisms and Associated UV-C Kill Dosage (99%)

PATHOGEN	BIOSAFETY LEVEL	UV Dose 99% (μW-s/cm ²)	PATHOGEN	BIOSAFETY LEVEL	UV Dose 99% (μW-s/cm ²)	PATHOGEN	BIOSAFETY LEVEL	UV Dose 99% (μW-s/cm ²)
Acinetobacter	2	3,600	Ebertelia typhosa	1	4,100	Proteus mirabilis	2	1,600
Adenovirus	2	11,800	Echovirus	2	1,600	Pseudomonas aeruginosa	1	10,500
Aeromonas	2	2,300	Eurotium (rubrum)	1	86,800	Reovirus	2	54,000
Aspergillus	2	19,200 - 896,000	Fusarium (solani)	1	62,600	Rhizopus	2	34,600 - 896,000
Bacillus anthracis	2	8,700	Haemophilus influenzae	2	7,700	Rhodotulula (spp.)	1	224,000
Bacillus magaterium sp. (spores)	1	5,200	Influenza A virus	2	6,600	Sarcina lutea	1	39,400
Bacillus magaterium sp. (veg)	1	2,500	Klebsiella pneumoniae	2	8,400	Scopulariopsis	2	578,000
Bacillus paratyphus	1	6,100	Legionella pneumophila	2	2,600	Serratia marcescens	1	21,000
Bacillus subtilis spores	2	11,000	Leptospira canicola - infectious Ja	1	6,000	Spirillum rubrum	1	8,800
Blastomyces dermatitidis	2	28,000	Listeria monocytogenes	2	31,100	Sporothrix schenckii	2	56,000
Botrytis cinerea	1	50,000	Measles virus	2	4,400	Staphylococcus albus	1	5,720
Burkholderia cenocepacia	1	11,600	Micrococcus candidus	1	12,300	Staphylococcus aureus	2	6,600
Candida albicans	1	150,000	Micrococcus sphaeroides	1	15,400	Staphylococcus epidermis	1	57,600
Cladosporium	2	37,800 - 896,000	Mucor (mucedo)	1	120,000	Staphylococcus hemolyticus	1	5,500
Clostridium perfringens	2	27,100	Mycobacterium avium	2	16,800	Staphylococcus lactis	1	8,800
Coronavirus	2	1,400	Mycobacterium kansasii	2	16,000	Streptococcus pyogenes	2	7,500
Corynebacterium diphtheriae	2	6,500	Mycoplasma pneumoniae	2	1,700	Streptococcus viridans	2	3,800
Coxsackievirus	2	23,000	Neisseria catarrhalis	2	8,500	Trichophyton	2	112,000
Cryptococcus neoformans	2	56,000	Nocardia asteroides	2	56,000	Ustilago (Zaeae)	1	224,000
Curvularia lunata	1	112,000	Phytomonas tumefaciens	1	8,500	Vaccinia virus	2	143,000
Molds			BIOLOGICAL AGENTS			Protozoa		
Aspergillus flavus	2	99,000	Hepatitis A	2	8,000	Chlorella Vulgaris	1	22,000
Aspergillus glaucus	2	88,000	Salmonella typhi	2	15,200	Paramecium	1	200,000
Aspergillus niger	2	330,000	Shigella	2	4,200			
Mucor racemosus A	2	35,200	Vibrio cholerae	2	6,500	Virus		
Mucor racemosus B	2	35,200				Bacteriophage - E. Coli	1	6,600
Oospora lactis	1	11,000	Yeast			Poliovirus - Poliomyelitis	2	6,000
Penicillium expansum	2	22,000	Brewers yeast	1	8,800	Tobacco mosaic	1	440,000
Penicillium roqueforti	2	26,400	Common yeast cake	1	13,200			
Penicillium digitatum	2	88,000	Saccharomyces cerevisiae	1	13,200			




Microgravity Science Glovebox

Microgravity Science Glovebox (MSG)



Marshall Space Flight Center

Decontamination System

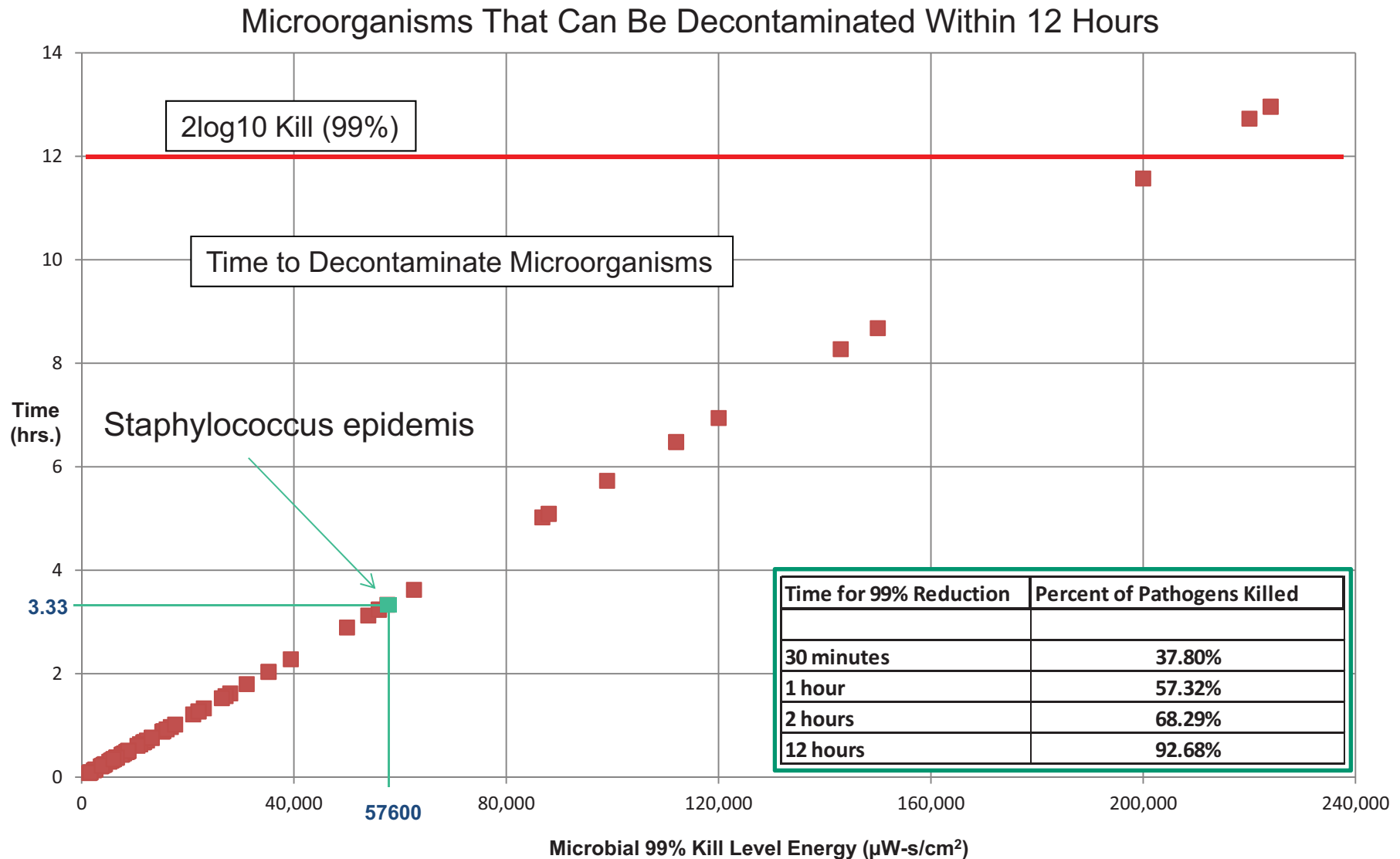



Microgravity Science Glovebox

Microgravity Science Glovebox (MSG)



Dexterous/Tactile Gloves

- **Biotech Gloves**
 - Thinner Gloves that provide more dexterity and sense of touch
 - 7 mil Hypalon Glove
 - Typical exam gloves are ~6 mils
- **Will adapt existing MSG design**



MSG has four glove ports; two on the front window and one on each side port. Glove ring assemblies can be installed in any glove ports as required by an investigation.

Gloves will be provided in three sizes
7,9, & 10.



MSG Glove & Gauntlet Configuration
(7 mil Hypalon Glove, 15 mil Gauntlet)



MSG Iris & Gauntlet Configuration



MSG Iris & Gauntlet With Nitrile Disposable Glove



Microgravity Science Glovebox

Microgravity Science Glovebox (MSG)



Marshall Space Flight Center

Video Upgrade Equipment (VUE) Available in 2014




Microgravity Science Glovebox

Microgravity Science Glovebox (MSG)




Marshall Space Flight Center

Current MSG Video System

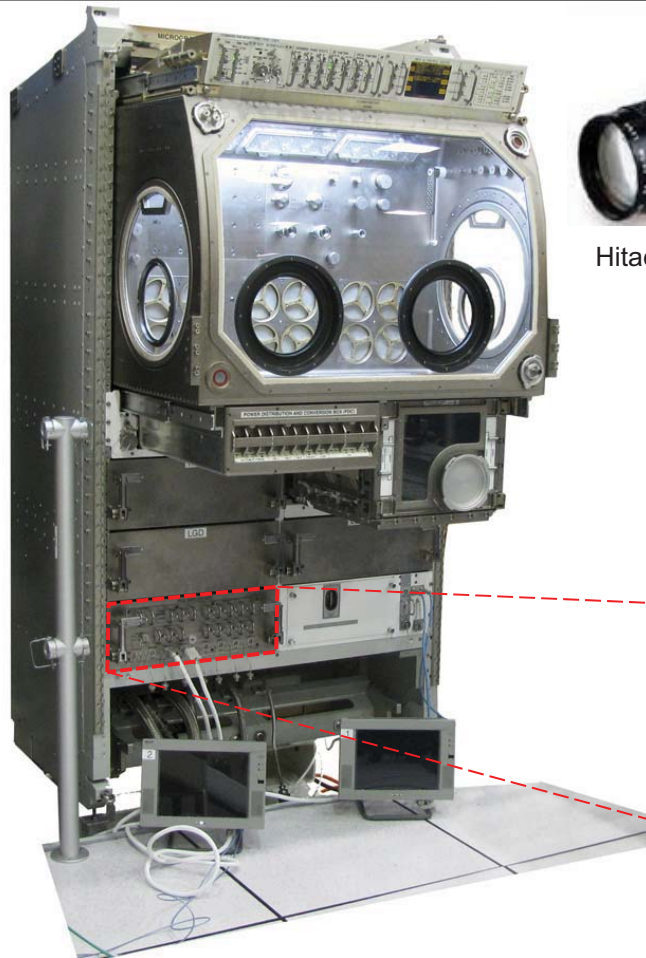

Hitachi HV-C20 Color Camera

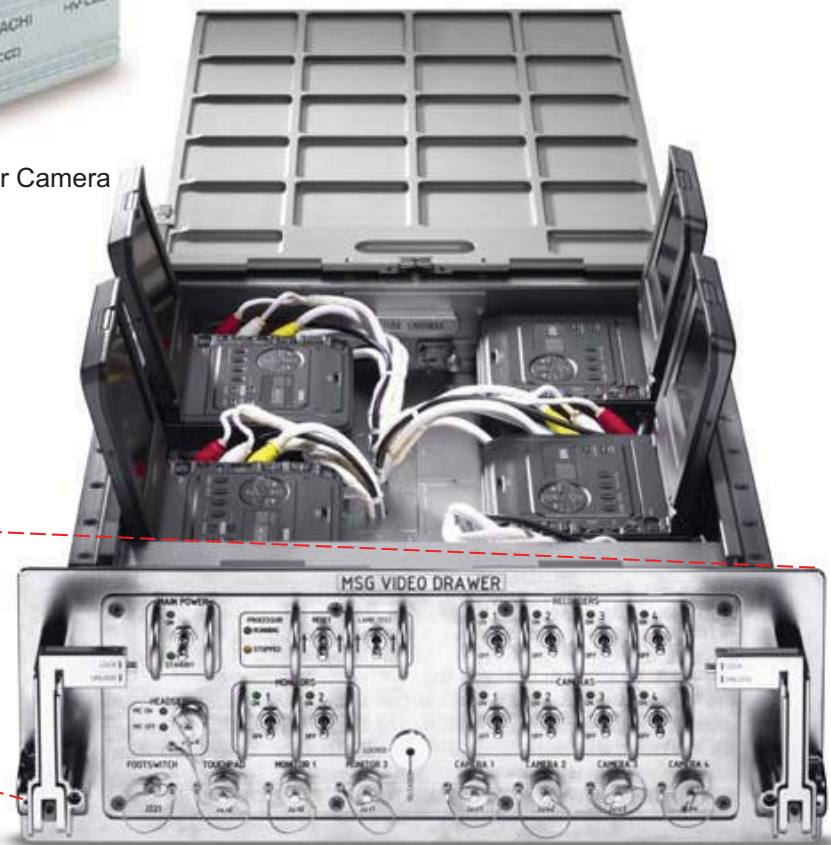

Pictured above in the bottom left drawer location of the MSG Engineering Unit, the MSG Video Drawer is shown connected to two video monitors. The Video Drawer is the main component of the MSG Video System.

In addition to accommodating 4 exchangeable video recorders, the Video Drawer contains power, communications, and remote control systems. The front panel allows for the crew to switch power to individual cameras, recorders, and monitors and to connect the various external components, including cameras and monitors.



Typical MSG Video System Setup





Microgravity Science Glovebox

Microgravity Science Glovebox (MSG)

Video System Overview



Marshall Space Flight Center

-
- The MSG Video Upgrade Equipment (VUE) will be capable of recording, storing, and transferring high definition/high resolution/high speed, color digital video data to ISS for downlinking.
 - The VUE will utilize significantly higher video resolution and speeds than the existing MSG video system thereby enhancing research observation activities
 - The MSG VUE consist of the following enhancements:
 - Powered ISIS drawer containing computer control and supporting electronics
 - High speed/high resolution cameras
 - High definition video cameras
 - GigE compatibility
 - Six terabytes of data storage via two 2 Tb Solid State RAID drives and two 1 Tb conventional hard drives.
 - Digital video data output capabilities for ISS to ground downlink. Downlink rates - up to 6 Mbps or higher depending on available bandwidth of the ISS LAN.



Microgravity Science Glovebox

Microgravity Science Glovebox (MSG)



VUE Camera Summary

Name	Type	Resolution	Sensor Size	Max Output
Prosilica 1050C	GigE	1024H x 1024V	1/2" Type CCD	1024 x 1024 w/ 8/12 Bit Color up to 109 fps
Prosilica 1910C	GigE	1920H x 1080V	2/3" Type CCD	1920 x 1080 w/ 8/12 Bit Color up to 55 fps
Flare 2KSDI	HD-SDI	2048H x 1088V (1920H x 1080V)	2/3" Type CMOS	2048 x 1088 w/ 10 Bit Color up to 30 fps
Hitachi HV C20 (Existing – to be replaced)	Analog RGB	768H x 494V	1/2" CCD	768 x 494 @30fps



Microgravity Science Glovebox

Microgravity Science Glovebox (MSG)



VUE Cameras

GX1050



GX1910



Size w/o lens (inches)
1.7 L x 2.5 W x 2.5 H
(w/o connectors)

Shown with Non-VUE Lenses



Microgravity Science Glovebox

Microgravity Science Glovebox (MSG)

VUE Hardware Description Cameras



Camera Info:

- Flight configuration: Two HD-SDI (Flare) cameras & Two Gig-E (Prosilica) cameras
- Two types of Gig-E cameras
 - 1910C 1920Hx1080V @ 56 fps
 - 1050C 1024Hx1024V @ 110 fps
- Each camera has a fixed, 10' long cable w/modified rear housing
 - HD-SDI camera will require a new feed through connector
 - This camera's cable is two headed

Note: Lenses are not installed on the depicted cameras



Prosilica

Flare



Microgravity Science Glovebox

Microgravity Science Glovebox (MSG) VUE Hardware Description Monitors



Monitor Info:

- Flight configuration utilizes two ViewPoint monitors
- Each monitor has a fixed, 10' long cable
- Monitors are for use external to the MSG Working Volume
- The hardware is MOTS



- 12.1" Wide Screen
- Resolution (1280x800 WXGA)
- Viewing Angle from all sides is 88 degrees
- 12VDC @ ~ 20 Watts

* Flight Monitor connectors are located on the bottom right of the units (as viewed from the front).




Microgravity Science Glovebox

Microgravity Science Glovebox (MSG) VUE Hardware Description Drawers



Drawer Info:

- Flight configuration is a single powered ISIS drawer
- Power is sourced through the rear drawer power connector and through a new J01 Jumper Cable
- Drawer was GFE to MSG by Boeing

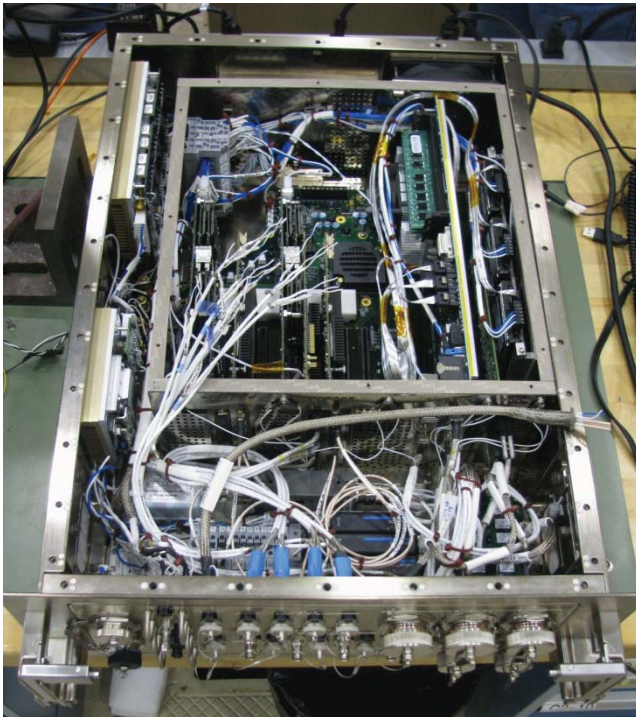

- Front panel interfaces include:
 - Power jumper and MLC
 - Cameras (8x) & monitors (2x)
 - Ethernet (3x)
 - USB (2x)
- Drawer is a standard 4 panel unit height
- Drawer & CPU tops are affixed w/threaded fasteners



Microgravity Science Glovebox



Conclusion

- The MSG is a very versatile and capable research facility on the ISS.
- The Microgravity Science Glovebox (MSG) on the International Space Station (ISS) has been used for a large body of research in material science, heat transfer, crystal growth, life sciences, smoke detection, combustion, plant growth, human health, and technology demonstration.
- MSG is an ideal platform for gravity-dependent phenomena related research. Moreover, the MSG provides engineers and scientists a platform for research in an environment similar to the one that spacecraft and crew members will actually experience during space travel and exploration.
- The MSG facility is ideally suited to provide quick, relatively inexpensive access to space for Physical Science, Life Science, and Biological Science investigations.





Microgravity Science Glovebox

Microgravity Science Glovebox (MSG)



Acknowledgements

The authors would like to thank the members of the MSG Integration Team for their assistance, and review of this paper. Specifically, Andrew Tygielski, Phillip Bryant, Chris Butler, Jeff Smith, Sharon Manley, and Mark Shelton provided invaluable information and assistance.